

STORZ MEDICAL



Before & After Pictures Acoustic Wave Therapy (AWT®)

Indication

Body shaping with Radial AWT® & diet program

Parameter per area

V-ACTOR® – 2500 pulses, 2.8 – 3.4 bar, 34 Hz

D-ACTOR® – 4000 pulses, 2.6 – 4.0 bar, 16 Hz, DI15

V-ACTOR® – 1800 pulses, 3.4 bar, 31 Hz

Before AWT®



Treatments

20

Clinic

Health & Vitality Salon, Netherlands

After AWT®



Indication

Body shaping with Radial AWT® & diet program

Parameter per area

V-ACTOR® – 2500 pulses, 2.8 – 3.4 bar, 34 Hz

D-ACTOR® – 4000 pulses, 2.6 – 4.0 bar, 16 Hz DI15

V-ACTOR® – 1800 pulses, 3.4 bar, 31 Hz

Before AWT®



Treatments

20

Clinic

Health & Vitality Salon, Netherlands

After AWT®



Indication

Body shaping with Radial AWT® & diet program

Parameter per area

V-ACTOR® – 2500 pulses, 3.0 – 3.8 bar, 32 Hz

D-ACTOR® – 4000 pulses, 2.6 – 4.0 bar, 16 Hz, DI15, D-20S

V-ACTOR® – 1800 pulses, 3.8 bar, 32 Hz

Before AWT®



Treatments

20

Clinic

Health & Vitality Salon, Netherlands

After AWT®



Indication

Body shaping with Radial AWT® & diet program

Parameter per area

V-ACTOR® – 2000 pulses, 3.0 – 3.6 bar, 32 Hz

D-ACTOR® – 3000 pulses, 2.6 – 3.4 bar, 16 Hz, C15

V-ACTOR® – 2000 pulses, 3.6 bar, 32 Hz

Before AWT®



Treatments

20

Clinic

Health & Vitality Salon, Netherlands

After AWT®



Indication

Body Shaping

Parameter per area

Radial AWT® – 4000 pulses, 2.6 bar, Deep Impact

Radial AWT® – 4000 pulses, 3.4 bar, D20-S

Treatments

10

Clinic

Cosmetic Centre, Great Britain

Before AWT®



After AWT®



Indication

Body Shaping

Parameter per area

Radial AWT® – 4000 pulses, 4.2 bar, Deep Impact DI15

Treatments

10

Clinic

Fitness Studio GALAXY, Austria

Before AWT®



After AWT®



Indication

Body Shaping

Parameter per area

Radial AWT® – 4000 pulses, 4.2 bar, Deep Impact DI15

Treatments

10

Clinic

Fitness Studio GALAXY, Austria

Before AWT®



After AWT®



Indication

Body Shaping

Parameter per area

Planar AWT® – 1000 pulses, 0.38mJ/mm²

Radial AWT® – 2000 pulses, 3.0 bar, D20-S

Treatments

8

Clinic

Dr. Goisis, Italy

Before AWT®



After AWT®



Indication

Body Shaping & Contouring

Parameter per area

Planar AWT® – 2000 pulses, 0.76mJ/mm²

Radial AWT® – 3000 pulses, 3.0 bar, D20-S

Treatments

6

Clinic

Dermacenter Stockholm, Sweden

Before AWT®



After AWT®



Indication

Body Shaping

Parameter per area

Planar AWT® – 1000 pulses, 0.76mJ/mm²

Radial AWT® – 2000 pulses, 1.8 bar, D20-S

Treatments

6

Clinic

Dr. Rümmelein, Switzerland

Before AWT®



After AWT®



Indication

Body Shaping & Cellulite

Parameter per area

Planar AWT® – 1000 pulses, 1.24 mJ/mm²

Radial AWT® – 4000 pulses, 3.0 bar, D15

Treatments

10

Clinic

Dr. Ledermann, Switzerland

Before AWT®



After AWT®



Indication

Body Shaping & Circumference Reduction

Parameter per area

Planar AWT® – 1000 pulses, 0.38mJ/mm²

Radial AWT® – 2000 pulses, 3.0 bar, D20-S


Treatments

11

Clinic

ANU Therapy, Ireland

Before AWT®

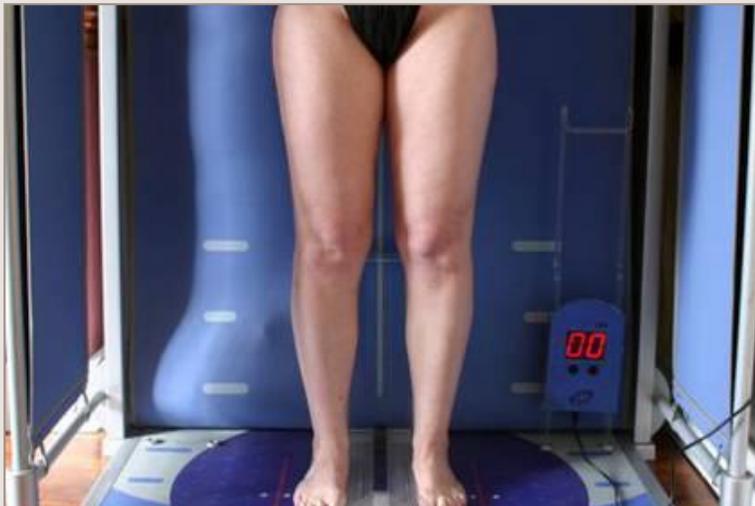

After AWT®



Indication

Cellulite

Parameter per area

Radial AWT® – 3000 pulses, 2.6 bar, D20-S

Treatments

8

Clinic

Dr. Johnson, USA

Before AWT®



After AWT®



Indication

Cellulite

Parameter per area

Radial AWT® – 3000 pulses, 2.6 bar, D20-S


Treatments

8

Clinic

Dr. Johnson, USA

Before AWT®

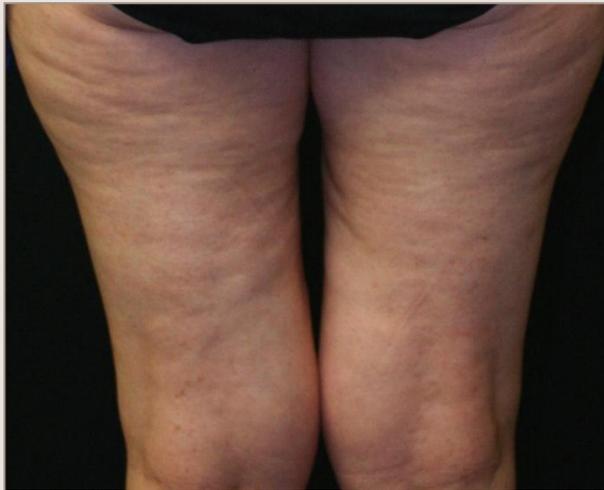

After AWT®



Indication

Cellulite

Parameter per area

Radial AWT® – 3000 pulses, 2.6 bar, D20-S

Treatments

8

Clinic

Dr. Johnson, USA

Before AWT®



After AWT®



Indication

Body Shaping & Cellulite

Parameter per area

Planar AWT® – 1000 pulses, 0.76mJ/mm²

Radial AWT® – 2000 pulses, 1.8 bar, D20-S

Treatments

6

Clinic

Dr. Rümmelein, Switzerland

Before AWT®



After AWT®



Indication

Cellulite

Parameter per area

Radial AWT® – 3000 pulses, 2.6 bar, D20-S

Treatments

8

Clinic

Dr. Johnson, USA

Before AWT®



After AWT®



Indication

Post Liposuction

Parameter per area

Planar AWT® – 2500 pulses, 0.76mJ/mm²

Radial AWT® – 4000 pulses, 2.2-3.0 bar, DI15

Radial AWT® – 1000 pulses, 2.2-3.0 bar, D20-S

Before AWT®



Treatments

6

Clinic

Dermacenter, Sweden

After AWT®



Indication

Post Liposuction

Parameter per area

Planar AWT® – 2500 pulses, 0.76mJ/mm²

Radial AWT® – 4000 pulses, 2.2-3.0 bar, DI15

Radial AWT® – 1000 pulses, 2.2-3.0 bar, D20-S

Before AWT®



Treatments

10

Clinic

Dermacenter, Sweden

After AWT®



Indication

Scar smoothing

Parameter per area

Planar AWT® – 3000 pulses, 0.76mJ/mm²

Radial AWT® – 3000 pulses, 3.0 bar, D20-S

Treatments

6

Clinic

Bodyform, Denmark

Before AWT®



After AWT®



Indication

Pregnancy stretchmarks, Striae

Parameter per area

Radial AWT® – 1000 pulses, 1.8 bar, Deep Impact DI15

Radial AWT® – 2000 pulses, 2.6 bar, D20-S

Treatments

6

Clinic

Cellulite-Praxis, Germany

Before AWT®



After AWT®

